


HandCut

Manual tru-cut needle for soft tissue biopsy with guillotine system

The HandCut-Needle is a disposable device for soft and fibrous tissues biopsy that need to be cut for their removal. The “Menghini” sharpening of the cannula allows cutting and acquisition of specimen. The stylet has a triple ultra-sharp tip to penetrate easily into fibrous tissues as well and has a 20mm specimen notch for sample capturing.

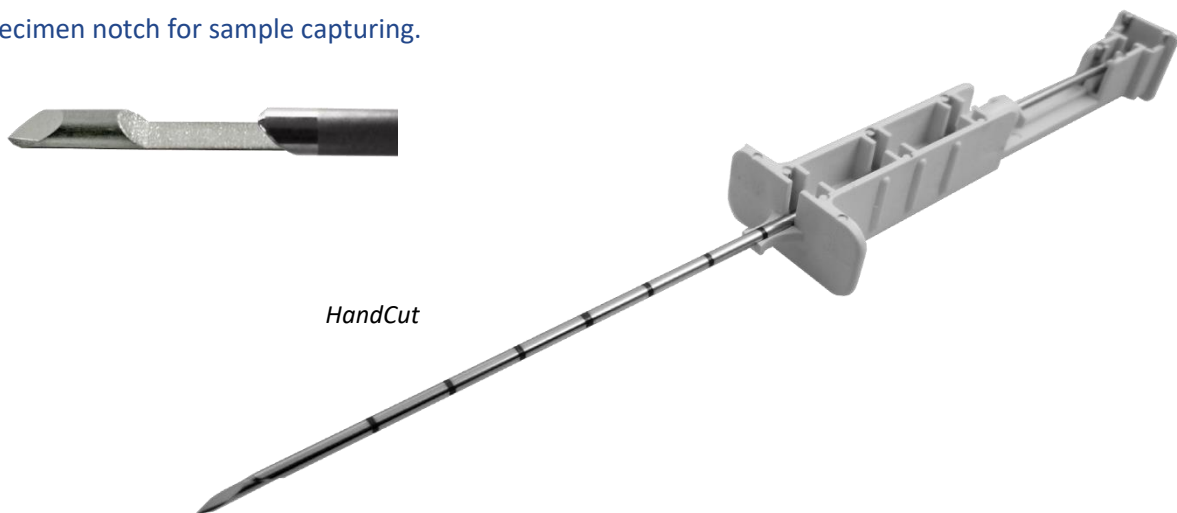

HandCut

- The device is equipped with a **echogenic marker** to enhance the ultrasound visibility avoiding cells spreading during the biopsy procedure
- The needle includes the **luer lock technology** that guarantees a secure connection to syringe for cytological aspiration
- The device is has a **length limiter** to take samples of 10 mm in length

Ordering Information:

Code	Gauge	Diameter mm	Length
PM14114	14 G	2.10	114 mm
PM1415	14 G	2.10	150 mm
Packaging unit: 10 pieces			

For other sizes, please contact our product specialist